

Econo-Cap & Heliflex GC Columns

Heliflex AT[®]-210+

- 他社相当品: DB-210, Rtx[®]-200, OV[™]-202, MEGA-200, OV[™]-210, OV[™]-215, QF-1, SP-2401, RSL-400

Heliflex AT[®]-210+ Specifications

Phase:	Trifluoropropyl Methylpolysiloxane
Polarity:	Polar
USP Code:	G6
Ideal For:	ケトン、アルデヒド、シラン、グリコール、ニトロ芳香族、除草剤、EPAメソッド 8140 & 609

Heliflex AT[®]-210+ Capillary Columns

Length	i.d.	Film	Temp. Limits* min.– max.	Old Reference	Part No.
15m	0.53mm	1.00µm	40– 240°C	13127	8635189
	0.53mm	1.20µm	40– 240°C	985115	8635190
30m	0.53mm	1.00µm	40– 240°C	13136	8635293
	0.53mm	1.20µm	40– 240°C	985130	8635294

Heliflex AT[®]-1701+

- 他社相当品: DB-1701, OV[™]-1701, SPB[™]-1701, Rtx[®]-1701, 007-1701, PAS-1701, CP-Sil 19CB, MEGA-1701, HP-1701, BP-10, ZB-1701

Heliflex AT[®]-1701+ Specifications

Phase:	(14% Cyanopropylphenyl)-86% Methylpolysiloxane
Polarity:	Intermediate Polarity
USP Code:	n/a
Ideal For:	殺虫剤、PCB's、ドラッグ、除草剤、TMS 糖類


Heliflex AT[®]-1701+ Capillary Columns

Length	i.d.	Film	Temp. Limits* min.– max.	Old Reference	Part No.
15m	0.25mm	0.25µm	-20– 280°C	13881	M8635566
	0.32mm	0.25µm	-20– 280°C	13762	M8635567
30m	0.25mm	0.25µm	-20– 280°C	13686	M8635558
	0.25mm	1.00µm	-20– 260/280°C	13687	M8635568
	0.32mm	0.25µm	-20– 280°C	13690	5121863
	0.32mm	1.00µm	-20– 260/280°C	13691	M8635517
	0.53mm	1.00µm	-20– 260/280°C	13795	M8635569
60m	0.53mm	1.20µm	-20– 260/280°C	13830	8635M405
	0.25mm	0.25µm	-20– 280°C	13688	M8635570
	0.25mm	1.00µm	-20– 260/280°C	13689	M8635571
	0.32mm	0.25µm	-20– 280°C	13692	M8635572
	0.32mm	1.00µm	-20– 260/280°C	13693	M8635573
	0.53mm	1.00µm	-20– 260/280°C	14131	M8635574

*: Maximum Temperature range may change depending on stationary phase film thickness.
In this catalog are reported only the most common GC columns dimensions and part numbers.
Contact us for other items and custom products.

EPA Method 552.1: Haloacetic Acids

1. Monochloroacetic Acid
2. Monobromoacetic Acid
3. Dalapon
4. Dichloroacetic Acid
5. Trichloroacetic Acid
6. Bromochloroacetic Acid
7. Dibromoacetic Acid

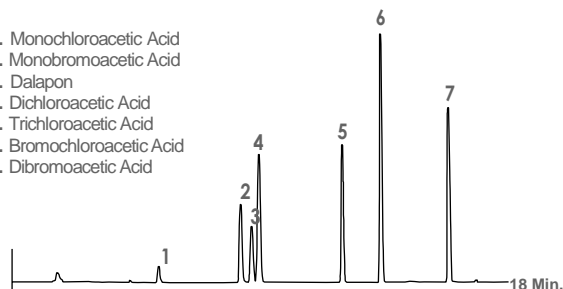

Column: Heliflex AT[®]-1701+, 30m x 0.32mm x 0.25µm
Part No. **5121863** (Old. Ref. 13690)
Temp: 50°C (10min hold) to 120°C at 10°C/min
Carrier Gas: Helium at 35cm/sec
Detector: ECD

Contact us for more info

Custom Capillary Coiling **Available**

Custom Capillary Column Dimensions **Available**

Custom Capillary Column with Integral Guard **Available**

FAST-GC narrow-bore Column Dimensions **Available**